

THE HUMAN BODY HAS FINALLY MET ITS MATCH


CUSTOM MOLDED SEATING SYSTEMS—ONLY FROM PRM.



SIGNATURE-FIT™



YOUR CLIENTS DON'T COME OFF AN ASSEMBLY LINE.

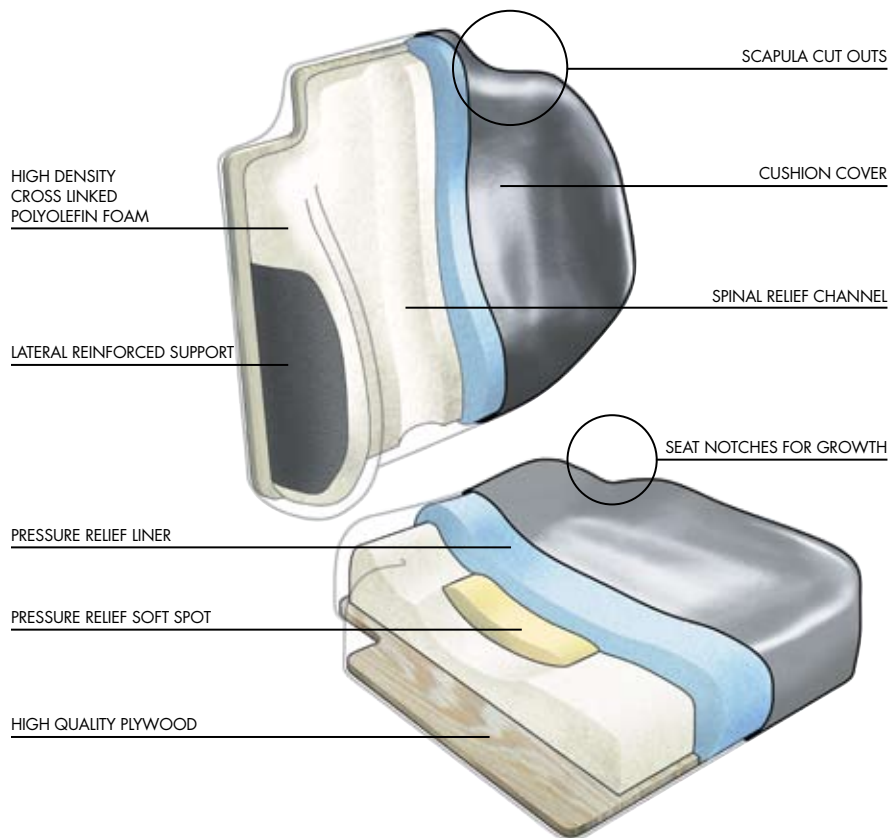

NEITHER DO OUR CUSHIONS.

AS YOU ARE WELL AWARE, NO TWO PEOPLE ARE EXACTLY ALIKE. THEY DIFFER IN PERSONALITY, LOOKS, HAIR COLOR, BODY TYPE AND SIZE. JUST TO NAME A FEW. WHICH IS WHY WE ARE PROUD TO PROVIDE OUR SIGNATURE-FIT CUSTOM MOLDED SEATING SYSTEM.



PRESSURE RELIEF LINER OPTIONS

All pressure relief liners are composed of high-density visco-elastic foam. The client shape is dilated to accommodate the thickness and firmness of the chosen liner before the carving process begins. This allows for compression of the liner and ensures proper fit. The most common foam chosen is soft in 1/2" or 1" thickness. We also offer medium firmness (1/2" and 1" thickness) and super soft liner foam (1/2" or 1"). Foam combinations can also be used.

CUSHION CONSTRUCTION

We mount the carved, high-density cross-linked polyolefin foam to Baltic birch plywood. The plywood has the T-nut pattern specified at time of order to facilitate ease of wheelchair attachment.

REINFORCED SOLID LATERAL SUPPORTS

Signature-Fit cushions provide extremely solid lateral support. This is accomplished by reinforcing the cushions using 1/8" or 1/4" ABS plastic. The plastic is attached to the side of the cushion and comes out to match the lateral depth—affording a solid lateral support while maintaining a thin profile. Additionally, a heavy-duty L-shaped reinforcement may be requested.

THE SIGNATURE-FIT DIFFERENCE.

ONCE YOUR CLIENT'S SHAPE IS DIGITIZED, YOU CAN SIMPLY E-MAIL THE DATA TO US FOR IMMEDIATE FABRICATION. SAVING BOTH VALUABLE TIME AND MONEY. DIGITAL SIGNATURE IS EASY TO USE, EXTREMELY ACCURATE AND A VALUABLE TIME SAVER FOR YOUR REHAB TECHNOLOGY SUPPLIERS. FOR ADDED VERSATILITY, WE CAN ALSO DIGITIZE DIRECTLY FROM CASTS OR FOAM MOLDS.



CUSTOM CUSHIONS

Signature-Fit custom cushions utilize both firm and soft foam to provide the added support your clients require—and the soft feel they demand. This foam combination maximizes the pressure-relieving qualities of the cushion shape.

CUSHION COVER OPTIONS:

We can finish cushions using any of the following materials. The best material is chosen depending on the environment and client needs.

VACUUM-FORMED FABRICS



VINYL

EASY TO MAINTAIN AND CAN BE WIPED CLEAN WITH SOAP AND WATER. OVER EIGHT COLOR CHOICES AVAILABLE.

DARTEX®

EASILY MAINTAINED FABRIC OFFERS STRETCH TO MAXIMIZE PRESSURE-RELIEF CHARACTERISTICS OF FOAM.



NEOPRENE®

PROVIDES FOUR-WAY STRETCH MAXIMIZING PRESSURE-RELIEF CHARACTERISTICS OF FOAM. AVAILABLE IN SIX COLORS.




SEWN REMOVABLE COVERS

AVAILABLE WITH OR WITHOUT A ZIPPER IN INCONTINENT-PROOF NEOPRENE® POLAR TECH® AND LYCRA® FOR CLIENTS WITH SIGNIFICANT INCONTINENCE ISSUES WE RECOMMEND TWO REMOVABLE COVERS.

EACH SIGNATURE-FIT CUSHION IS AS UNIQUE AS ITS OWNER.

TRIAL FIT POLICY

The Signature-Fit carving process enables shipment of unfinished cushions for trial fit. This provides you with the opportunity to make adjustments to the shape prior to final shipment. If no modifications are needed, the trial fit process can be bypassed and the cushions will be shipped within 7-10 working days.

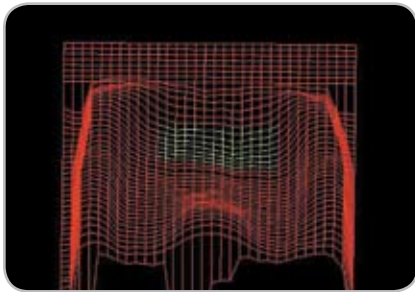

CUSHION MODIFICATIONS

Signature-Fit cushions can be fabricated to enhance the mold shape and/or facilitate attachment to any wheelchair. These options include pressure relief soft spots, rail-cuts, footrest removal notches, rear seat notches and more. Please see the order form for additional options.

BACK MOUNTING

Back mounting options include adjustable drop hardware. You may also specify which hardware you are using, and we will provide the applicable T-nuts for attachment.

SEAT MOUNTING

Seat mounting options include Velcro® and adjustable drop hardware. You may also specify which hardware you are using, and we will provide the applicable T-nuts for attachment.

SHORT LEAD TIMES

At PRM, we'll ship both trial fit and finished Signature-Fit cushions within 7-10 working days of order receipt. For more details, samples or to schedule a demonstration with a PRM representative, call 866-PRM-REHAB or click prmrehab.com today.

PRM. THE POWER OF PRECISION.

Your clients rely on you to provide them with comfortable, supportive seating that fits—and most importantly—that lasts. At PRM, our custom cushion selections do that and more. Precision-Fit offers maximum pressure relief, increased sitting tolerance, postural support and comfort finely tuned to each individual's positioning needs. You won't get that from just any off-the-shelf cushion. For details, samples or to schedule a demonstration with a PRM representative, call **866-PRM-REHAB** or click prmrehab.com today.

11861 East Main Road
North East, PA 16428
ph: 814/725-8731
fx: 814/725-2934
toll free: 866-PRM-REHAB
www.prmrehab.com



SIGNATURE-FIT™